

Certificate of Analysis

Veyronalabs
Veyrona Labs
Slovakia
support@veyronalabs.com

Liquilabs s.r.o.
Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Semax 10 mg
Batch Number CPT-SEMAX10-20260725-007
Date Published 2026-09-15 13:24

Results for LYO-0851

Peptides	Result	Unit	Uncertainty	Acceptable Range
Semax Assay Peptide Screening 0.1% TFA	11.33	mg	[± 0.06]	
Semax Purity Peptide Screening 0.1% TFA	98.8	%	[± 0.5]	
Semax Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	983		[± 5]	
Semax Identification by RT Peptide Screening 0.1% TFA	0.992		[± 0.005]	
Microbiology	Result	Unit	Uncertainty	Acceptable Range
Bacterial Endotoxin Chromgenic USP<85>/ Eur. Ph. 2.6.14. Bacterial Endotoxin Chromgenic Test	0.120	EU/mg	[± 0.002]	

	Method Specification	
Determination of identity, content and purity of Semax		
<i>Document number</i> SMX_008_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% DFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% DFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.8 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	225 nm


Gradient Program		
Time [min]	A [%]	B [%]
1	98	2
13.99	35	65
14	5	95
15	5	95
15.01	98	2
19	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

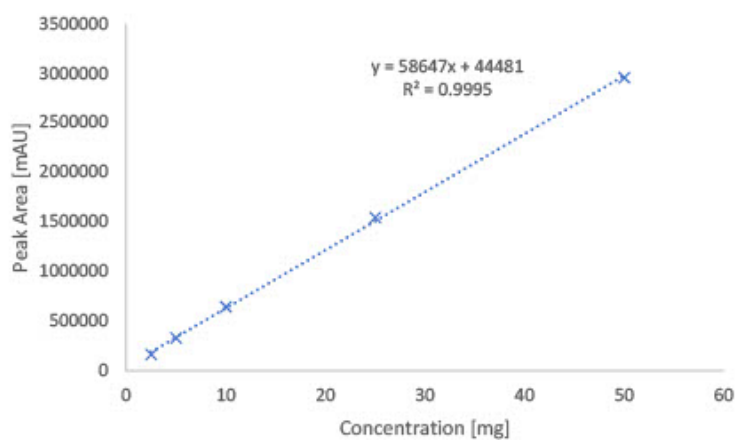

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% DFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% DFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.8 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Waters XSelect CSH C18, 100x2.1mm 2.5µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
1	98	2
13.99	35	65
14	5	95
15	5	95
15.01	98	2
19	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225 nm.

Analysis Report




Analysis of quantity and purity of active ingredient by UHPLC with UV detection

Sample Information


Injection Volume : 2
Data File : LYO-0851_012.lcd
Method File : Peptide Screening_DFA_Group E.lcm
Date Acquired : 9/10/2026 7:35:52 PM

Chromatogram

uAU

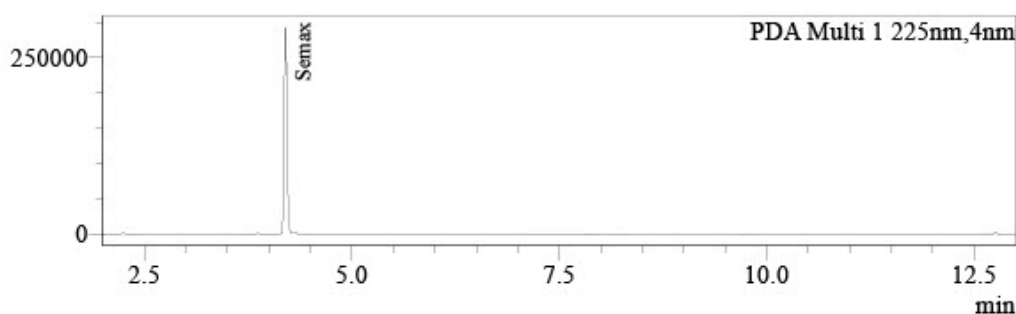

uAU

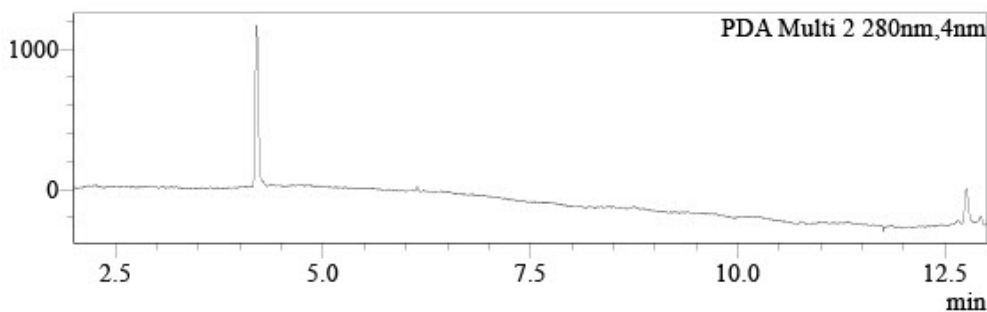

Peak Table

PDA Ch1 225nm

Name	Ret. Time	Area	Conc.	Unit	Area%
	3.675	411	0.000		0.057
	3.866	1398	0.000		0.195
	4.086	713	0.000		0.099
Semax	4.197	708891	11.329	mg	98.779
	4.318	684	0.000		0.095
	5.310	108	0.000		0.015
	5.663	221	0.000		0.031
	6.135	397	0.000		0.055
	12.759	4828	0.000		0.673
		717651			100.000

Peak Table

PDA Ch2 280nm

Name	Ret. Time	Area	Conc.	Unit
	4.200	3149	0.000	
	14.221	4356	0.000	
		7505		

Endotoxin Determination Report




File Information

Filename: C:\UVVis-Data\Data\File_260910_110406.vqud
Date/Time: 09/10/2026 11:59:39 AM

Instrument Information

Instrument Type: UV-1900 Series
Model (S/N): UV-1900i (A12536253123)

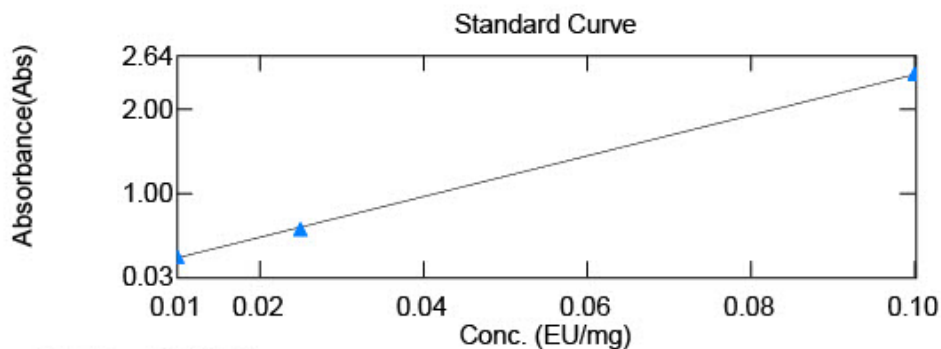

$$y = 24.1316x + 0.00350132$$
$$r^2 = 0.99988$$

	Sample Name	Conc	Raw_WL545.0	Dilution Rate
1	LYO-0838	0.119	0.1471	20.00
2	LYO-0839	2.461	1.4880	40.00
3	LYO-0840	0.120	0.1481	20.00
4	LYO-0841	0.123	0.1516	20.00
5	LYO-0842	0.145	0.1787	20.00
6	LYO-0843	0.119	0.1469	20.00
7	LYO-0844	0.634	0.7689	20.00
8	LYO-0845	2.332	1.4106	40.00
9	LYO-0846	0.076	0.0947	20.00
10	LYO-0847	1.766	1.0689	40.00
11	LYO-0848	0.062	0.0779	20.00
12	LYO-0849	0.081	0.1009	20.00
13	LYO-0850	2.660	1.6081	40.00
14	LYO-0851	1.122	1.3567	20.00
15	LYO-0852	0.081	0.1015	20.00
16	LYO-0853	0.657	0.7958	20.00
17	LYO-0854	0.637	0.7720	20.00
18	LYO-0855	2.207	1.3349	40.00
19	LYO-0856	0.076	0.0951	20.00
20	LYO-0857	0.080	0.1004	20.00

	Sample Name	Conc	Raw_WL545.0	Dilution Rate
21	LYO-0858	0.068	0.0853	20.00
22	LYO-0859	0.072	0.0906	20.00
23	LYO-0860	0.081	0.1008	20.00
24	LYO-0861	2.603	1.5736	40.00
25	LYO-0862	0.064	0.0813	20.00
26	LYO-0863	0.079	0.0985	20.00
27	LYO-0864	0.070	0.0877	20.00
28	LYO-0865	0.362	0.0908	100.00
29	LIQ-0068	0.007	0.1667	1.00

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.